


Customizing Vector

OpenFreezer uses the restriction sites designated on the insert as ligation points in the vector. The user can override this automatic creation by customizing the cloning sites. The user can also select to reverse complement the insert prior to insertion into the vector. Hybrid Sites can be created by selecting different but compatible enzymes for the vector and insert.

Customize Sites

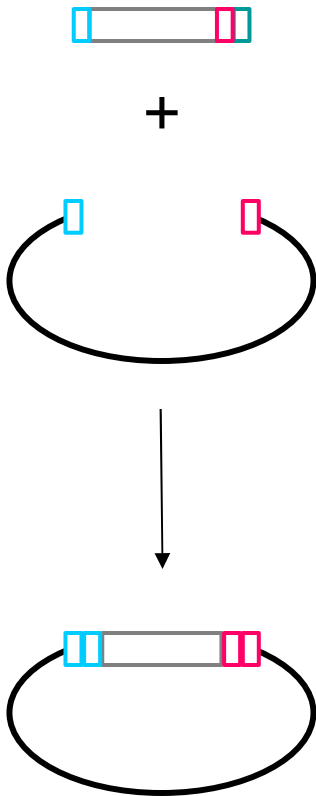


Customize Sites

5' Vector and Insert – BglII
3' Vector and Insert – EcoRI

(override BamHI as 3' ligation site)

Reverse Complement

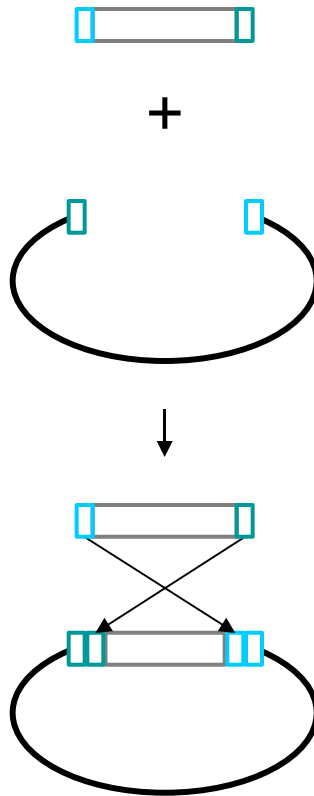


Customize Sites

5' Vector and Insert – BamHI
3' Vector and Insert – BglII

(designated 3' BamHI site on Insert is switched to a 5' BamHI site on new Vector)

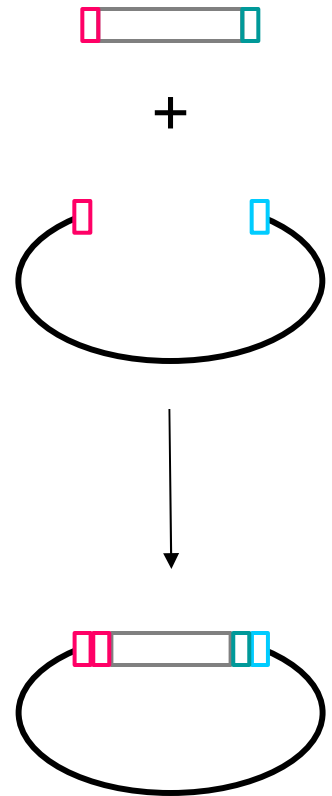


Create Hybrid Site(s)



Customize Sites

5' Vector and Insert - EcoRI
3' Vector - BglII, 3' Insert – BamHI

(OpenFreezer creates a Hybrid site at the 3' ligation site)

 BamHI
 BglII
 EcoRI